

Patient Name: _____ Contact Name: _____ Height: _____ Weight: _____ Date: _____

Ship Name & Address: _____ Account Name: _____ Account #: _____ P.O #: _____

Product Information


Product includes one Legassist Calf and One Compreboot Plus Foot.

Right Leg Foam Additional Foot Option
 Left Leg Flat Wavefoam Custom Medaboot (additional charge)

Note: If the greatest circumference measurement is >60cm, order a Super.

Circumference Measurements

(Lateral aspect of leg)

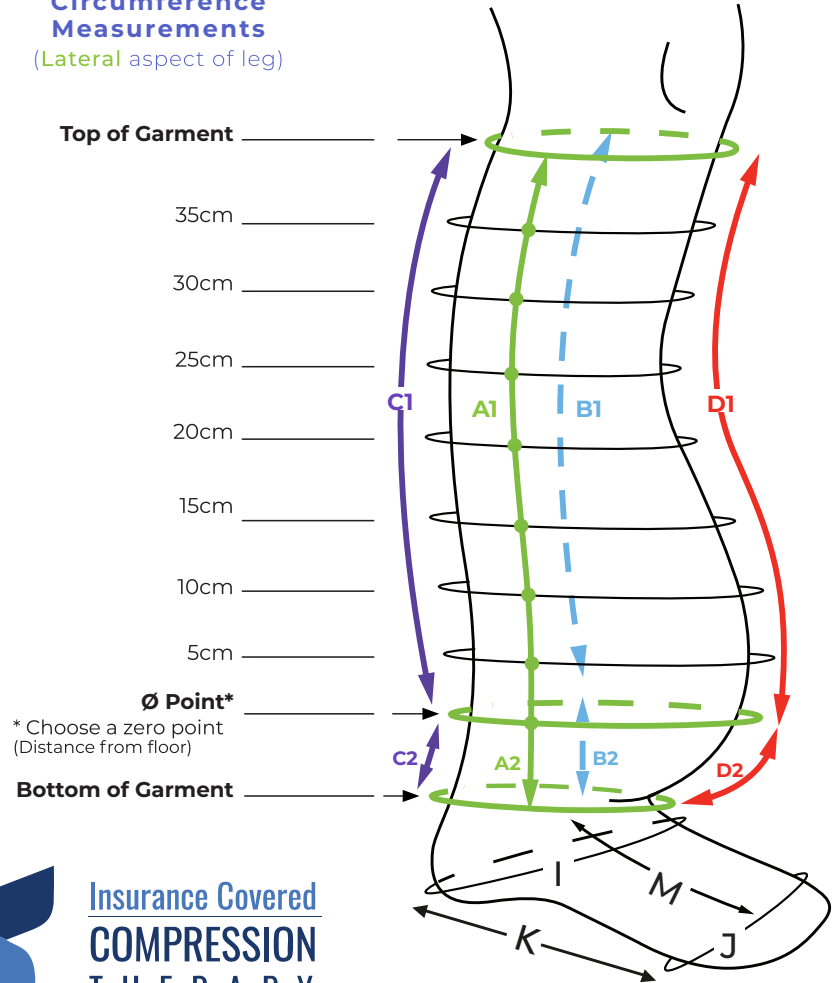

Length Measurements

Lateral

A1 _____
A2 _____

Medial

B1 _____
B2 _____

Posterior

C1 _____
C2 _____

Anterior

D1 _____
D2 _____

Foot Measurements

I _____
J _____
K _____
M _____



Insurance Covered
**COMPRESSION
THERAPY**

Phone 309-664-7930 • Fax 309-664-7931
Email orders@medsourceorders.com

Scan the QR code
to view a Measuring
Instruction video



A PDF of this order form can be found
online at: sigvaris.group/mceforms
or scan here:



CUSTOM Measuring Instructions

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our MCE Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera. photos of the leg with measurement markings must be emailed to: **us_orders@sigvaris.com**
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

INSTRUCTIONS

All measurements should be recorded in centimeters. Apply slight tension to hold the tape measure in place.

STEP 1: Choose The Zero Point (Ø point)

The zero point is unique to the Legassist Calf & Foot With Shelf garment. It is the reference point from which all other measurements are taken, in order to determine the lengths needed for the custom garment to fit the patient's individual ankle shape. Encircle leg with bottom edge of tape measure immediately above the shelf, at the lowest aspect of the curve of the ankle lobule where a tape measure can be applied without sliding off. Trace a line around the leg along the bottom edge of the tape using the SIGVARIS GROUP body pen, and record the circumference on Ø point line. Document your zero point location by measuring the distance from the floor.

STEP 2: Choose Where The Garment Will End Distally - "bottom of garment"

Encircle leg with bottom edge of tape measure where the ankle lobule touches the ankle and/or foot. Do not pull tape under shelf. Mark leg on this edge of tape at the anterior, medial, posterior and lateral aspects, and record the circumference on the line labeled bottom of garment. The smallest ankle measurement is not relevant for this product. The product is designed to hug and lift the ankle lobule.

STEP 3: Choose Where The Garment Will End Proximally - "top of garment"

The ideal placement for the top edge of the garment is two finger widths below the popliteal crease, and across the fibular head and tibial tuberosity. However, the position of the top of garment can be customized, if needed, to overlap better with a Legassist LCS thigh piece. Encircle leg with top edge of tape measure at the desired location. Using the SIGVARIS GROUP body pen, trace a line around the leg along the top edge of the tape and measure and record the circumference on top of garment line.

STEP 4: Measure And Record Lengths

Measure each of the lengths below from the Zero Point, following the contour of the leg and ankle, but not extending the tape measure into folds. Lengths measured distal to the Zero Point should end where the ankle lobule touches the ankle or foot at the bottom of garment versus extending the tape measure up under the lobule. Record lengths as follows:

Lateral (A) Lengths:

A1: Measure & record lateral length from Ø point to top of garment

A2: Measure & record lateral length from Ø point to bottom of garment

Medial (B) Lengths

B1: Measure & record medial length from Ø point to top of garment

B2: Measure & record medial length from Ø point to bottom of garment

Posterior (C) Lengths

C1: Measure & record posterior length from Ø point to top of garment

C2: Measure & record posterior length from Ø point to bottom of garment

Anterior (D) Lengths

D1: Measure & record anterior length from Ø point to top of garment

D2: Measure & record anterior length from Ø point to bottom of garment

A PDF of this order form can be found online at: sigvaris.group/mceforms or scan here:



CUSTOM Measuring Instructions**STEP 5: Measure And Record Circumferences**

On the lateral aspect of the leg only, mark every 5cm moving proximally from the zero point. Next, record the circumferences at each point, making sure that the measuring tape is parallel to the floor.

STEP 6: Measure Foot

1. Measure the circumference from where the shelf touches instep of the foot to the apex of the heel and record on line I.
2. Measure the foot circumference at the metatarsal heads and record on line J.
3. Measure the length from the 1st metatarsal head to the posterior aspect of the heel, and record on line K.
4. Measure the length from the dorsal aspect of the third metatarsal head to where the shelf touches the instep of the foot and record on line M.

STEP 7: Pictures

Now that you have completed your measurements, take pictures of the leg (front, back, side) that show all of the measurement markings, as well as the overall shape of the leg. This will assist our expert design team to further customize the garment according to your patient's individual needs. E-mail the pictures with your completed forms to: us_orders@sigvaris.com. Don't forget to record the patient's age, height, and weight, as this also helps with customizing the garment.

STEP 8: Select Additional Options

You are ready to complete your forms and place the order.

Choose right or left leg

Choose Size: If the greatest leg circumference is >60cm, order a Super.

Choose Flat or Wavefoam (recommended)

Choose Compreboot Plus (included) or Custom Medaboot (additional charge)

TAKE DIGITAL PHOTOS (REQUIRED) AFTER MEASUREMENT MARKINGS ARE MADE. EMAIL TO: us_orders@sigvaris.com

NOTE: Legassist Calf & Foot With Shelf is available in:

7099-BK Flat Foam, 7199-BK WaveFoam, 7299-BK Super, Flat Foam, 7399-BK Super, WaveFoam.

A PDF of this order form can be found
online at: sigvaris.group/mceforms
or scan here:

